

SFE Analysis & Non-destructive contact angle measurements

Portable Contact Angle Meter PCA-11

A compact and lightweight device for fully automated on-site measurements of contact angle, surface free energy of solids, and surface tension

Testing on-site for quality control and evaluation of surface wettability, characterized by the contact angle, is often required during manufacturing processes or in the field.




The PCA-11 can be used on any material, such as glass, metal, film, paper and pulp, polymer, and wood. It is ideal for on-site measurements or samples that cannot be analyzed with the usual stationary contact angle measuring instruments due to their size or shape.

Weighing only 650g, combined with the single-click fully automated operation and precise measurements, the PCA-11 is the ideal solution for quickly analyzing solids' surface free energy and for non-destructive surface evaluations.



FEATURES

- Single-click, fully automated operation
- Precise measurement of contact angles, even on slopes of up to 90 degrees
- Quick analysis of surface free energy by a simple exchange of filled syringes
- Measurement of surface tension using the pendant drop method
- Time-dependent measurements by image capture of 30 fps
- One prepared syringe allows for 500 measurements of droplets with a volume of 1 μ l
- Simple replacement of syringe and liquid sample via the easily accessible front door
- The visible needle tip allows precise capture of the target point

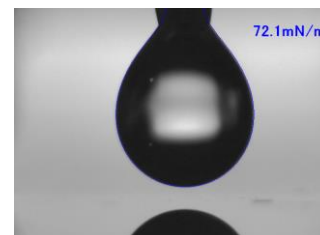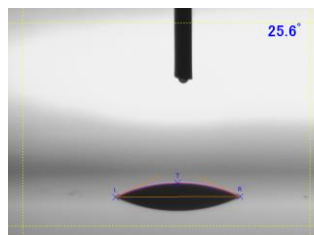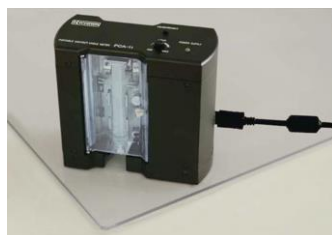

APPLICATIONS

- Evaluation of cleanliness of large panels and wafers
- Evaluation of surface modification techniques, such as corona, flame, and plasma treatments
- Evaluation of hydrophobic and hydrophilic surfaces from automobile parts, windshields, building walls, etc.

SPECIFICATIONS

Measuring methods	Contact angle: Sessile drop	Surface tension: Pendant drop
Analysis methods	Contact angle: Surface tension: Surface free energy of solids:	$\theta/2$, tangent, curve fitting (ellipse, circle), Young-Laplace Young-Laplace, ds/de OWRK, Owens-Wendt, Kaelble-Uy, Kitazaki-Hata, Wu, acid-base, Zisman
Measuring range	Contact angle: 0 to 180 °	Surface tension: 0 to 1000 mN/m
Resolution	Contact angle: 0.01 °	Surface tension: 0.01 mN/m
Accuracy ¹⁾	Contact angle: 0.5 °	Surface tension: 0.5 mN/m
Droplet dispensing	Automatic	
Dispensing resolution	0.1 μ l	
Droplet deposition	Automatic	
Power supply	USB bus power (a Windows laptop is required to operate the instrument)	
Dimensions, weight	110mm(W) x 50mm(D) x 110mm(H), 650 g	
Operating environment	Temperature: +10 to +35 °C, relative humidity: 30 to 80 % (non-condensing)	

¹⁾ Accuracy is the repeatability in terms of standard deviation based on the manufacturer's calibration standard.

For detailed information, please get in touch with our sales partner or us directly at +81-48-483-2629 or overseas-sales@face-kyowa.co.jp.

Specifications and designs are subject to change without notice.

2409